

A woman with long dark hair, wearing a grey athletic top and grey underwear, is shown from the waist up. The image is partially covered by a semi-transparent blue overlay on the left side. The background is a bright, outdoor setting, possibly a beach or a sunny area.

**CM**  
Slim  
Contouring Master

MAKE  
YOUR  
BODY  
SLIM

[www.CMSLIM.com](http://www.CMSLIM.com)

CE  
0068





CMSLIM is the most advanced and intensive electromagnetic muscle stimulator.

Focused electromagnetic field passes through all skin and fat in order to effectively stimulate the muscle, providing the most intensive continuous contractions which is ideal for muscle growth, and also inducing apoptosis.

CMSLIM is a HI-EMT device designed for aesthetic purpose, having 2 (two) applicators with higher intensity.

Its cutting-edge technology in non-invasive body contouring, as it not only burns fat, but also builds muscle.

This is key to a toned and healthy physique.

# New Technology in Non-invasive Body Slimming




Two applicators are placed on the target muscle area, such as your abs, thighs, or buttocks. The applicators then generate intense electromagnetic that causes involuntary muscle contractions.

# Two Applicators

These contractions trigger the release of free fatty acids, which break down fat deposits and increase muscular tone and strength.

This is very similar to the way doing workout.

As working out, the stimulation causes to rebuild and repair muscle tissue, resulting in stronger, firmer muscles.

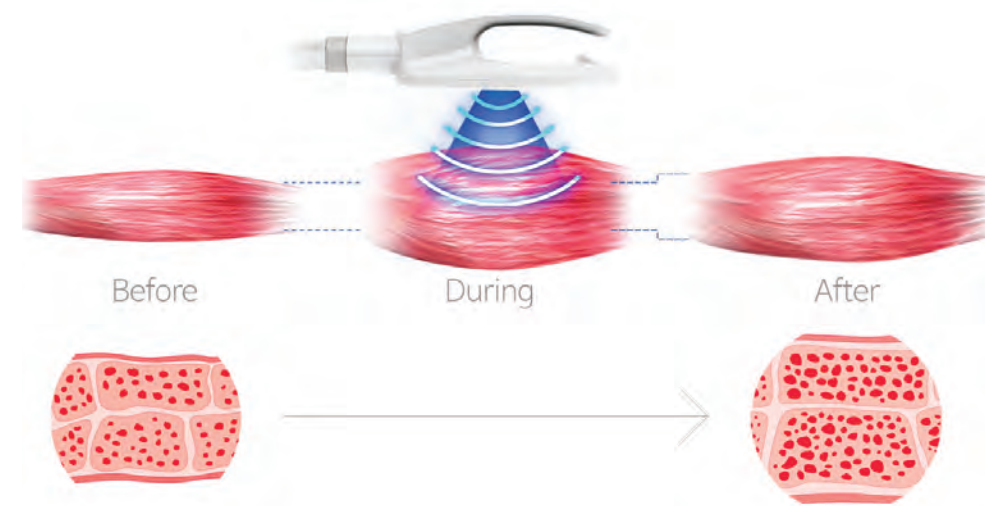




Non-invasive medical technology that is used for strengthening and re-education of muscles via interaction of the magnetic field with the tissue of the patient.


Contrary to voluntary muscle contractions, the supramaximal contractions are independent of brain function.

The HI-EMT uses a specific range of frequencies that does not allow muscle relaxation between two consecutive stimulations.

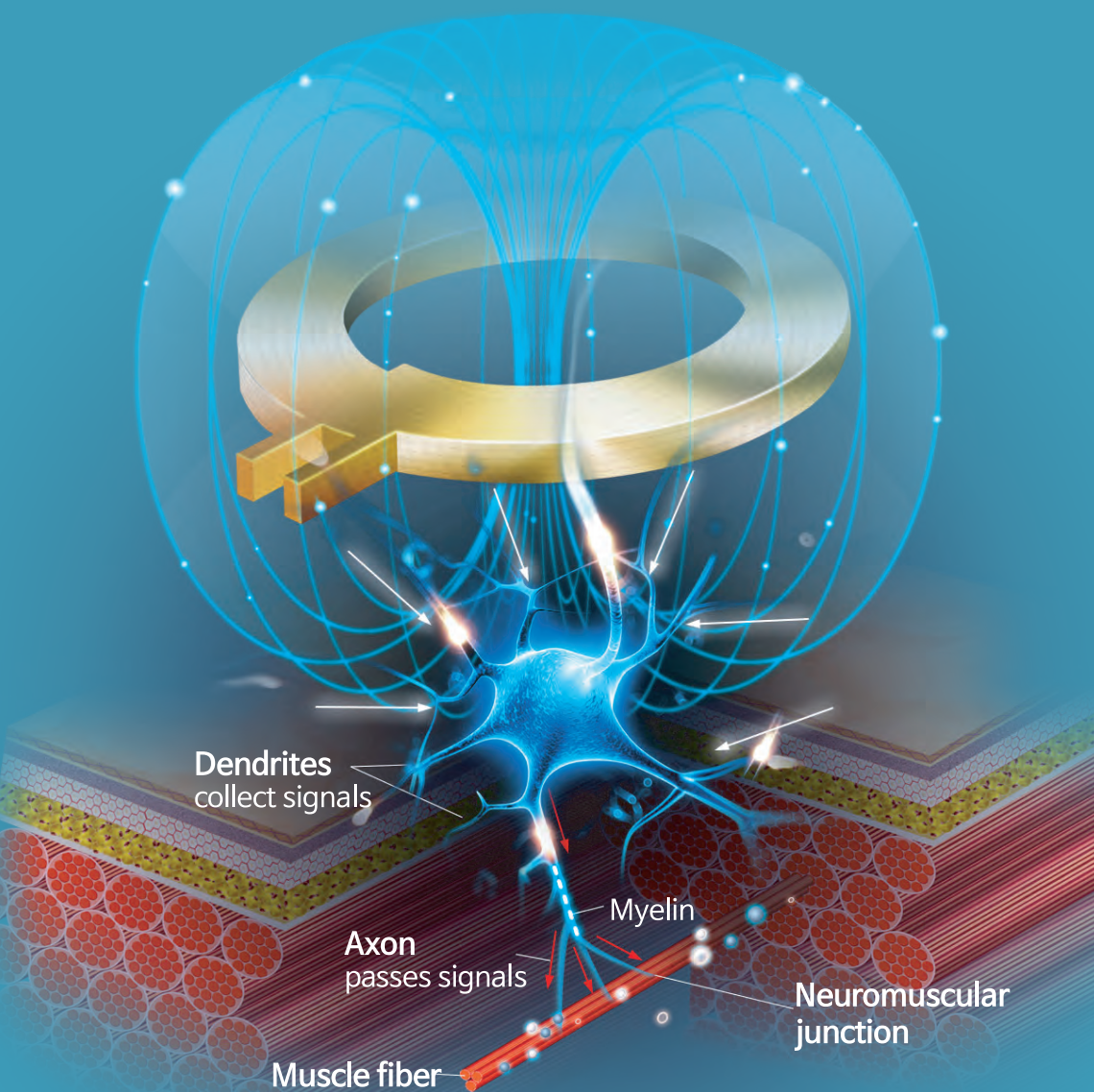

# Proven Technology

The HI-EMT uses a specific range of frequencies that does not allow muscle relaxation between two consecutive stimulations. The muscle is forced to remain in contracted state for multiple seconds. When repeatedly exposed to these high load conditions the muscle tissue is stressed and is forced to adapt.

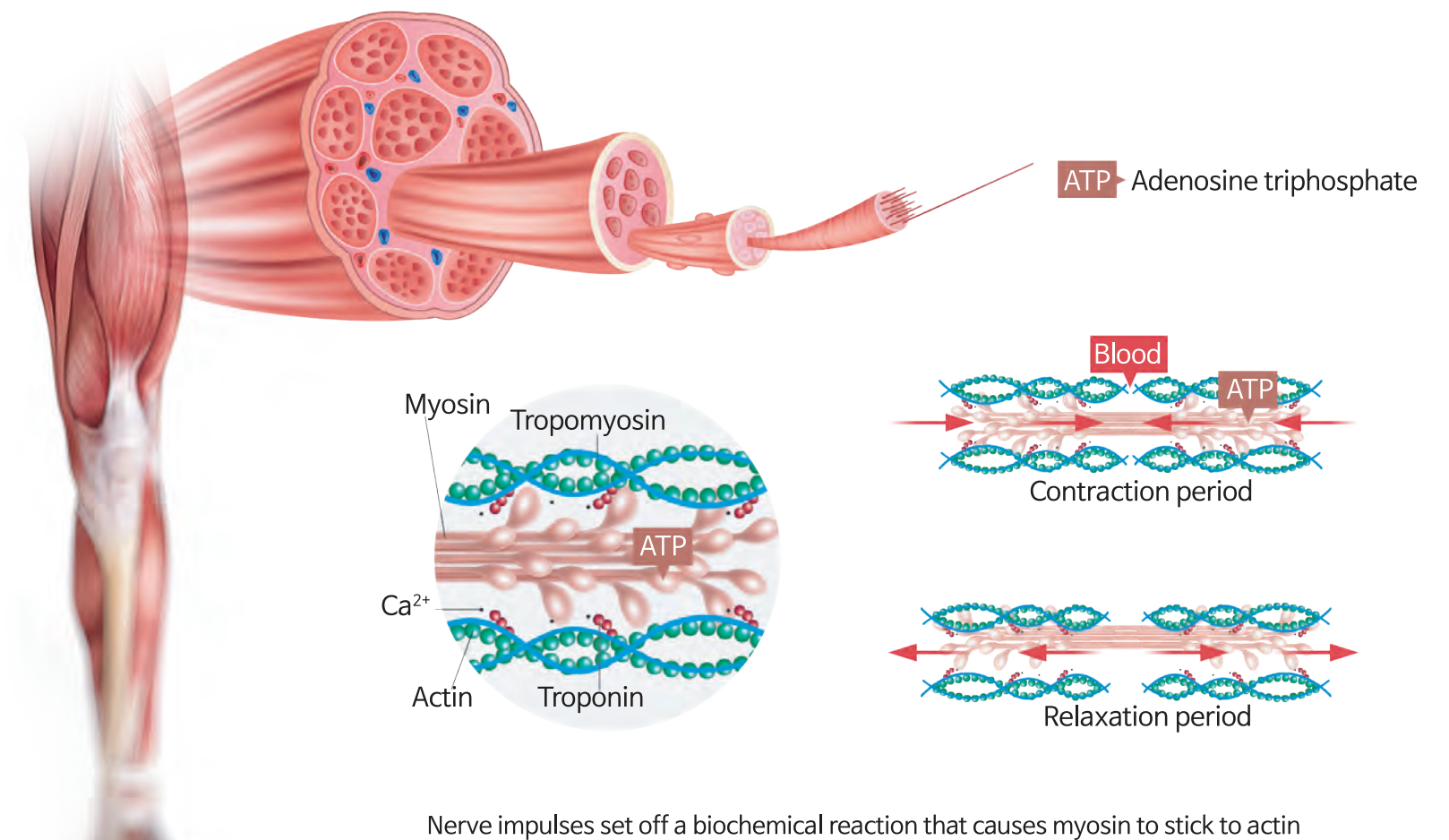

Mechanism of muscle contraction

Recent studies reported that on average 15% - 16% increase in abdominal muscle thickness was observed in treated patients one to two months after HI-EMT treatments.





# HI-EMT

## High Intensity Electromagnetic Muscle Trainer



is the medical technology that is used in aesthetic medicine. It uses focused electro- magnetic field with safe intensity level.

Electro- magnetic field passes non-invasively through the body and interacts with motor neurons which subsequently trigger supramaximal muscle contractions. “HI-EMT is used as a non-invasive body contouring treatment for abdomen and buttocks.”



Effect on Fat Several recent studies using CT, MRI and ultrasound evaluations have reported approximately 19% reduction in subcutaneous fat layer in patients treated by HI-EMT based device on their abdomen. The principle of cell apoptosis induced by increased concentrations of free fatty acids has been previously observed and demonstrated in numerous research studies.



# Session Information

CMSLIM's treatment session requires no anesthesia, incisions, or discomfort. In fact, patients are able to sit back and relax, while the device performs the equivalent of more than 20,000 painless crunches or squats.



**Hercules Course**  
Building Muscle

Treatment	1	2	3	4	5	6
Program	Hypertrophy	Combo2	Hypertrophy	Combo2	Hypertrophy	Combo2

**David Course**  
Burning Fat

Treatment	1	2	3	4	5	6
Program	HIIT	Combo1	Strength	Strength	HIIT	Combo1


Single treatment session takes 30 minutes only, and there is no downtime. 2-3 treatments per week, would be enough for perfect result for most people. Generally, 4-6 treatments recommended.





# Workout Program

CMSLIM's unique workout program consists of a series of an actual physical training program, which has been ideally designed for the maximum stimulation leading to the best result.

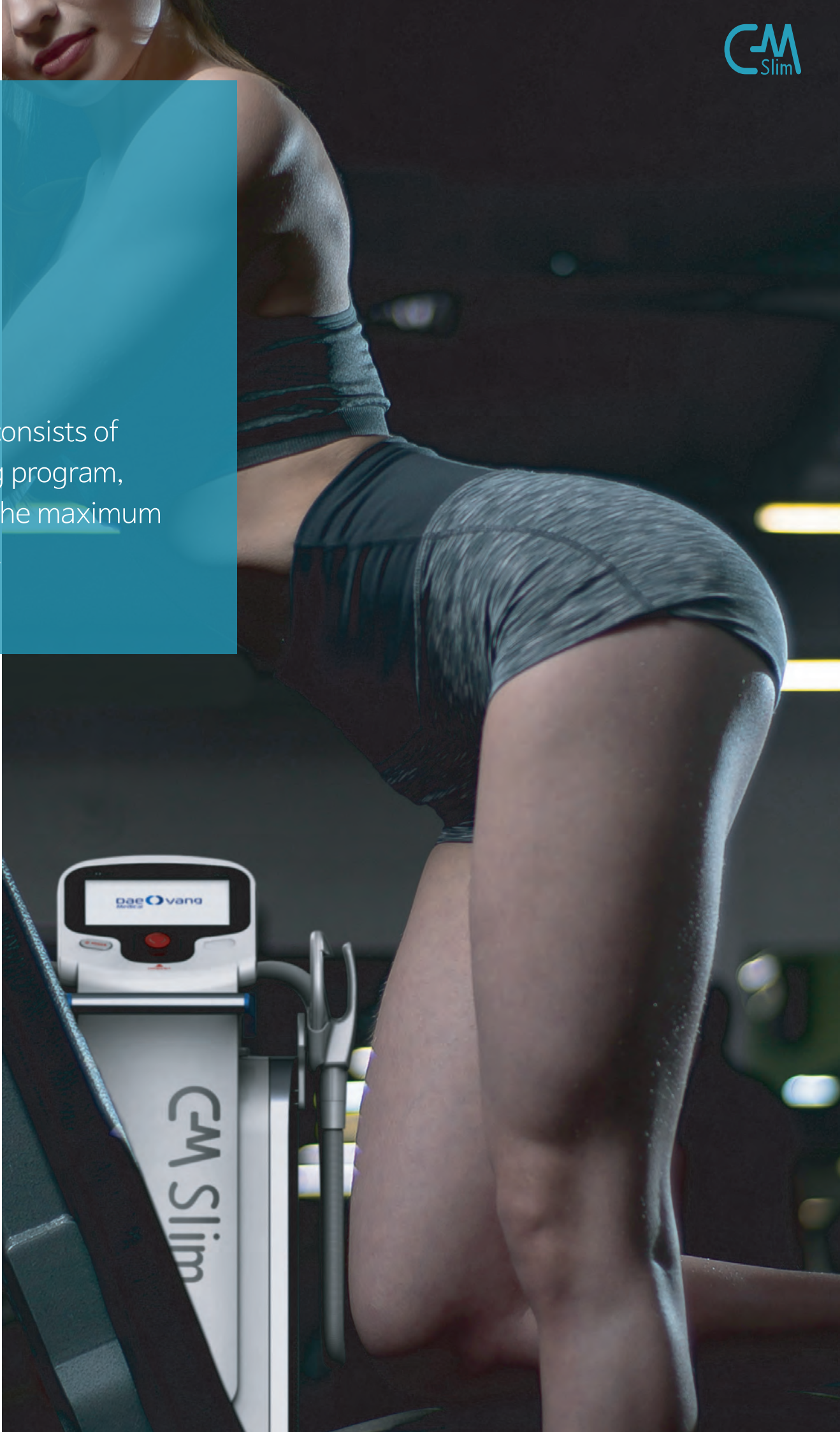

**Venus Course**  
Building Muscle

Treatment	1	2	3	4	5	6
Program	Hypertrophy	HIIT	Combo2	Hypertrophy	HIIT	Combo2


**Aphrodite Course**  
Burning Fat

Treatment	1	2	3	4	5	6
Program	HIIT	Combo 1	Strength	HIIT	Combo2	HIIT

Additionally, CMSLIM is an extremely efficient treatment that takes place over the course of 10-14 days, rather than months, and it also produces a more toned abs and uptight buttocks within two-to-four weeks.

If you need more information about our workout program, please visit [www.CMSLIM.com](http://www.CMSLIM.com) and let us know your inquiries.



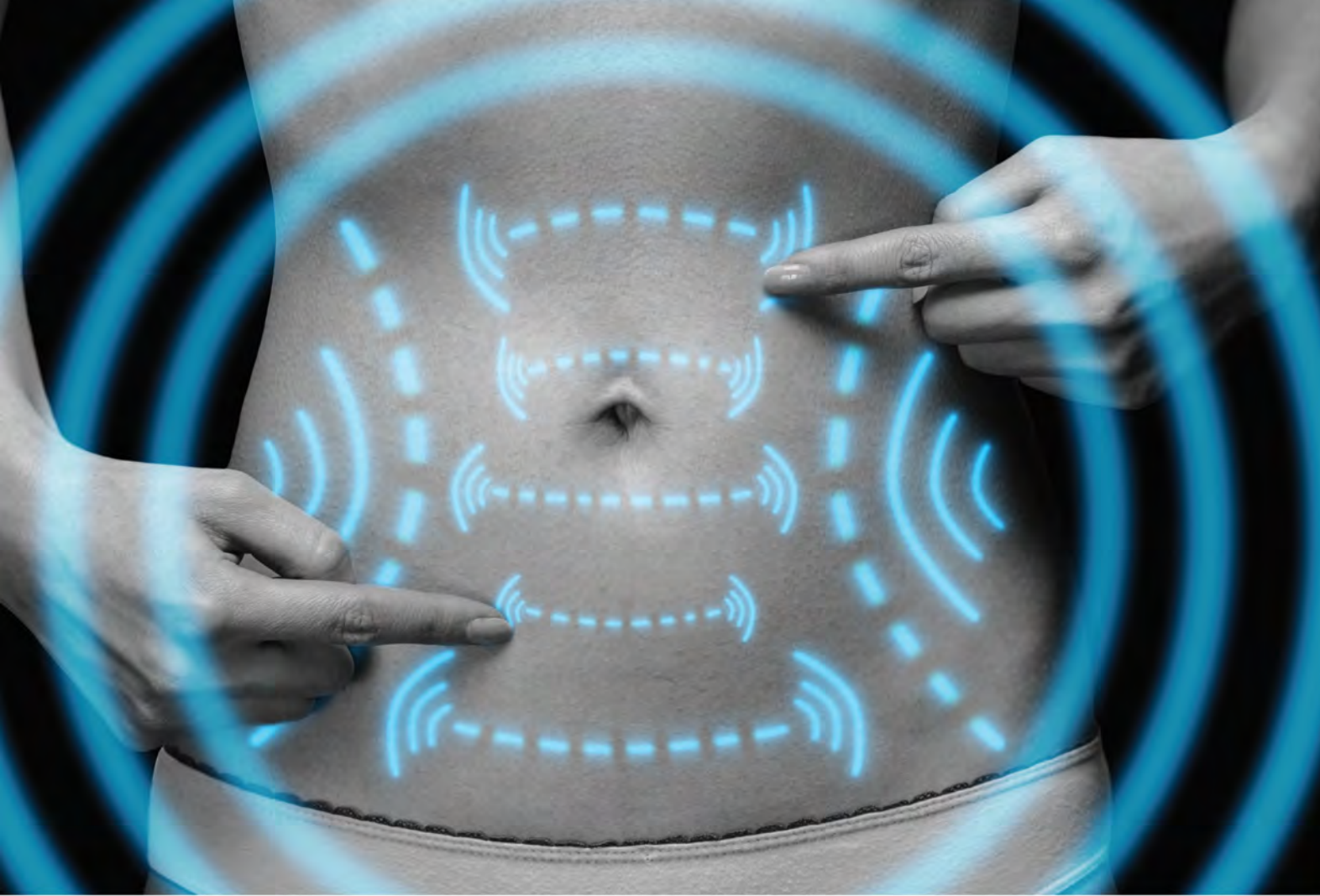

# CMSLIM Advantages


Output Intensity  
**2,5 Tesla**

HI-EMT base 2,5 Tesla output enables that electromagnetic energy could embrace big skeletal muscles of the human body. That means CMSLIM's muscle re-education process is deep enough like normal muscle contraction. Superficial electrical muscle stimulation is incomparable.



## Cooling Technology

Our optimized water cooling system make CMSLIM very effective to operate for a very long time without any downtime. Even further there are no worries for leaking oils anymore.

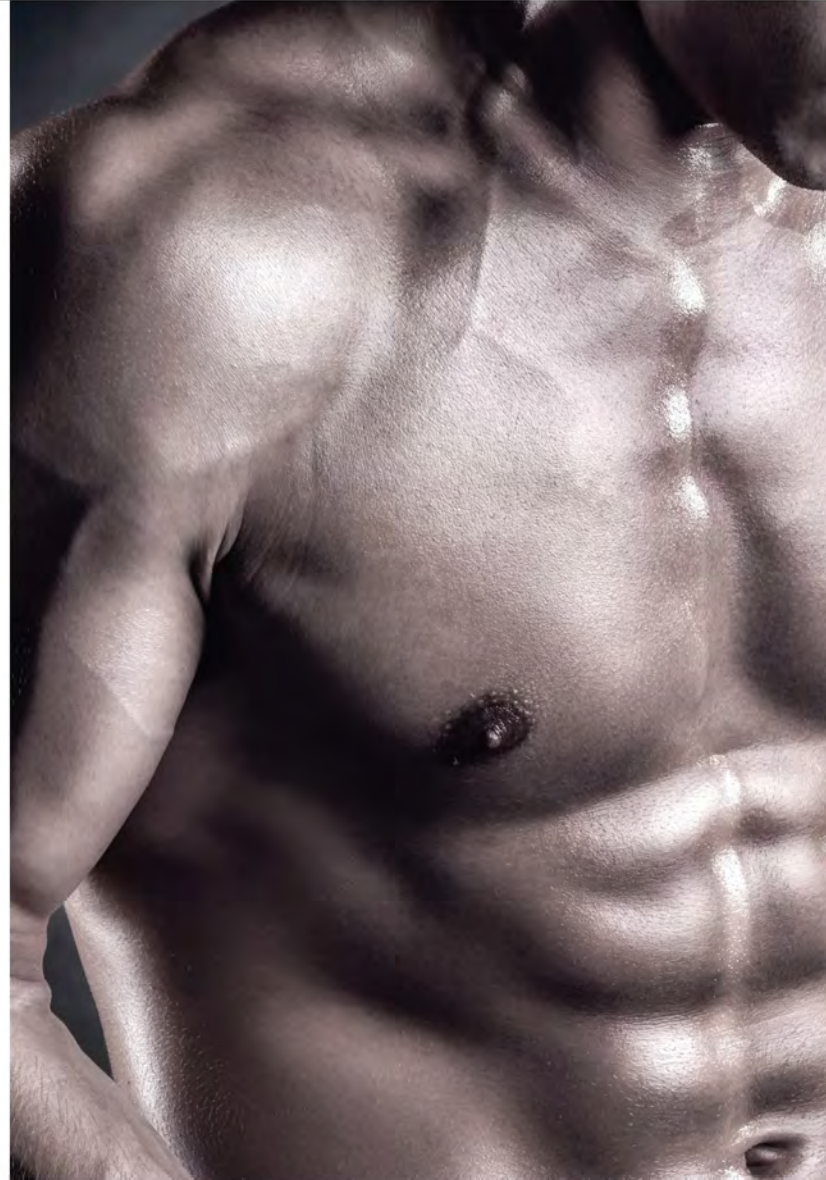

## Manual Training

Our cutting-edge technology makes it possible to adjust electromagnetic pulse width through very wide range. So you don't have to stick to a few presets, you can design special manual treatment program for your patient's needs.





No More Extra Supplies, CMSLIM has NO CONSUMABLES.



Perfect Water Cooling System makes CMSLIM EVEN STRONGER.



Great Durability – 10,000 Treatment Assurance



1 Year Full Coverage Warranty



Convenient Operation – EASY & EFFECTIVE PRESET / User Interface



# Maintenance



# Before & After

Even if you aren't in optimal shape, you can experience benefits from CMSLIM for years to come with regular exercise and balanced eating.



M1\_Before\_Front



M1\_After\_Front



M1\_Before\_Side



M1\_After\_Side



F1\_Before\_Front

F1\_After\_Front



F2\_Before\_Front

F2\_After\_Front



# FAQs

Q A candidate?

A **Anybody**

Q A downtime?

A **Not at all**

Q How long?

A **30 minutes per session**

Q How often?

A **2~3 times a week**

Q How many time?

A **4~6 treatment recommended**

Q How long can I use it?

A **No limit**

Q How long will it last?

A **Extremely long**

Q Additional supplies?

A **No consummables**

Q Warranty?

A **1 year**

Q Manual program?

A **Easy & Effective manual mode**



- ① LCD    ② Emergency stop button    ③ Reset Switch    ④ Applicator (2EA)

# Specification



Power	Power Input	230 V 50/60 Hz
	Maximum Power Consumption	2.3 kVA
Magnetic Stimulation	Magnetic induction amplitude (Intensity)	2,5 Tesla
	Shape of stimulation pulse	Biphase
	Pulse duration	300 $\mu$ s
Dimension	Size (height x width x depth)	420 x 660 x 1150 (mm)
	Weight	70 kg
	Safety type	BF
Country of Origin	Republic of KOREA	
Manufacturer	DAEYANG MEDICAL CO., LTD	



daeOyang  
Medical



[www.bioskin.it](http://www.bioskin.it)  
[shop.bioskin.it](http://shop.bioskin.it)



FACEBOOK



Follow Up

[www.dyat.life](http://www.dyat.life)  
[www.CMSLIM.com](http://www.CMSLIM.com)  
[www.daeyangmed.com](http://www.daeyangmed.com)

T. +82 (0) 33 747 4466 F. +82 (0) 33 747 1133  
147 Donghwagongdan-ro, Wonju-si,  
Gangwon-do, Republic of Korea  
[sales@daeyangmed.com](mailto:sales@daeyangmed.com)