


## Innovation of Dual Fractional Rejuvenation & Shining the Skin

### Total Fractional Capabilities

As the newest aesthetic trend, the fractional technology of ACTION II allows user-selected precise ablation depth and coagulation control with a painless procedure. The specialized multiple stacking technology of ACTION II supplies adjustable ablation depth and thermal coagulation zone.

The dynamic moving technique is available to shorten the operation time.

The fractional approach is available in both "Fractional" and "Long Pulse" modes, and controllable collateral thermal damage provides a virtually painless treatment. ACTION II fractional technology leads to a rapid healing process and minimal downtime.

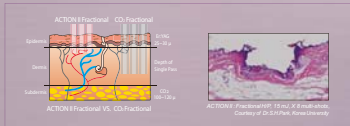

Specifications	
Laser Type	Er:YAG, 2940nm
Operating Mode	Normal Mode Long Pulse Mode Fractional Mode
Beam Profile	Top Hat
Pulse Energy	Normal Mode: 22J Fractional Mode: 18mJ
Pulse Width	Normal Mode: 250µs Long Pulse Mode: 1ms-1000ms
Spot Size	Zoom HWP: Normal: 1.2-7mm Ø Long Pulse: 1.2-5.4mm Ø Shooting Peel HWP: 14mm Ø Fractional HWP: 12mm Ø, 30mm Ø
Handpiece Auto Setting	Q
Pulse Rate	Single, 1, 2, 3 and 10Hz
User Interface	8 inch LCD Display / Touch screen
Memory	Shooting Peel: 80 spots
Aliming Beam	500mV-650mV (5400, 11000, Adjustable Brightness)
Dimension (mm)	250(W) x 650(D) x 1320(H)
Weight	86.1kg



Technology for Beauty and Health  
**LUTRONIC**

**Lutronic Corporation**  
Room 403-2, 2-45, 454 Nam Technobun 1141-1, baekok-dong  
Bundang-gu, Gyeonggi-do, Gyeonggi-do, Korea  
TEL: +82-31-908-3440 / FAX: +82-31-907-3440  
e-mail: [office@lutronic.com](mailto:office@lutronic.com)

**Lutronic Inc.**  
3003 N. 1st Street Suite 215 San Jose, CA 95134, USA  
TEL: +1-888-588-5644

**Lutronic Europe**  
122 Gillet Road, Edgbaston/Birmingham B16 9ES UK  
TEL: +44(0)121-454-3807 / Mobile: +44(0)7815487548

**Lutronic Japan Co., Ltd.**  
Mikusa Bldg., 3-21-47 Utsunomiya, Bunkyo-ku, Tokyo, 113-0034 Japan  
TEL: +81-3-3354-4227 / FAX: +81-3-3354-4228

**Lutronic China**  
Suite 22, Hyundai Motor Tower, 38 Kaoyun Road, Chaoyang District, Beijing, CHINA  
TEL: +86-10-84585940

Copyright©2011, Lutronic Corporation. All Rights reserved. Product and specifications are subject to change without notice. AC2B6PE1107

Technology for Beauty and Health  
**LUTRONIC**

## Easy & Optimal Multiplexing Platform

**ACTION II**  
4 in 1 Innovation of Er:YAG Aesthetic Laser



## Fractional Rejuvenation with Dual Operation Modes

## No More Limitations to Conventional Resurfacing and Coagulation

Designed for easily approach versatile aesthetic applications ACTION II is the cutting edge EY-AG laser platform of choice for a wide range of aesthetic procedures with the newly adopted Long Pulse Mode.

ACTION II indications include wrinkles, acne scars, large pores, skin tone, texture improvement and the removal of benign epidermal lesions, with additional possibilities as new treatment techniques are explored.

The latest ACTION II features the conventional micro-pulse mode for micro laser resurfacing and the long pulse mode for coagulation. The maximum 1000 ms pulse width in the long pulse mode is expanding new applications. Specialized fractional technology and the Shining Peel technique supply more gentle fractional rejuvenation and deliver a real glow to the skin.

ACTION II is equipped with Lutronic's proprietary Controlled Depth Ablation (CDA) technology and Top-Hat beam profile allowing precision control of ablation depth layer-by-layer with minimal collateral heat diffusion. This function maintains a consistent fluence for uniform results. This versatile laser system is good for both conventional micro laser resurfacing and deep resurfacing combined with the long pulse coagulation procedure. Properly controlled residual heating minimizes added discomfort, making procedures with ACTION II virtually painless.

For the fractional approach, Shining Peel and various indications, ACTION II is the ideal choice, reliably delivering controllable thermal energy with less patient discomfort, reduced treatment times and less patient downtime so that you may achieve higher patient satisfaction levels and improved revenue results.

## Shining Peel, the most gentle peeling technique for Fine Wrinkles, Large Pores and Skin Rejuvenation without downtime

### Shining Peel for Gentle Rejuvenation

Shining Peel is a brand new laser peeling approach for fine wrinkles, large pores and skin rejuvenation without topical anesthesia, downtime or and PIH. This gentle peeling technique stimulates collagen remodeling by delivering precise amounts of laser energy to the skin without immediate ablation of the epidermis. Shining Peel decreases the chance of PIH since the outside layer of the skin is kept intact during treatment. The outer epidermal layer will gradually flake off and this procedure helps patients avoid downtime or any other inconvenience after treatment. Shining Peel guarantees patients the shining face of youth.



**Shining Peel**  
Painless Treatment for Skin Rejuvenation

ACTION II is equipped with essential and convenient advantages such as the triple operation mode and the auto-detected triple handpieces for versatile clinical applications.

### Triple Operation Mode

- Normal mode
- Fractional mode
- Max 1000 ms Long Pulse mode - Fractional approach is available

### Adjustable Aiming Beam

The brightness of the 650 nm red diode aiming beam can be adjusted according to the lesion or skin color.

### Auto-Detected Triple Handpieces & Spot Size



1 - 7 mm Zoom Handpiece

### Smart User Interface Touch LCD

The smart user interface LCD touch screen and auto-detected handpiece provide convenient function.

## Major Indications of ACTION II

- Shining Peel
- Dual Fractional Rejuvenation : Large Pores, Skin Tone, Wrinkles
- Scars : Acne Scars, Surgical Scars
- Nevus-related & Epidermal Pigments  
Lentigo, Freckles, Seborrheic Keratosis (SK), Nevus, Syringomas, Warts

